

# Document made available under the Patent Cooperation Treaty (PCT)

International application number: PCT/US2004/030977

International filing date: 22 September 2004 (22.09.2004)

Document type: Certified copy of priority document


Document details: Country/Office: US  
Number: 60/504,774  
Filing date: 22 September 2003 (22.09.2003)

Date of receipt at the International Bureau: 02 May 2008 (02.05.2008)

Remark: Priority document submitted or transmitted to the International Bureau in compliance with Rule 17.1(a) or (b)



World Intellectual Property Organization (WIPO) - Geneva, Switzerland  
Organisation Mondiale de la Propriété Intellectuelle (OMPI) - Genève, Suisse

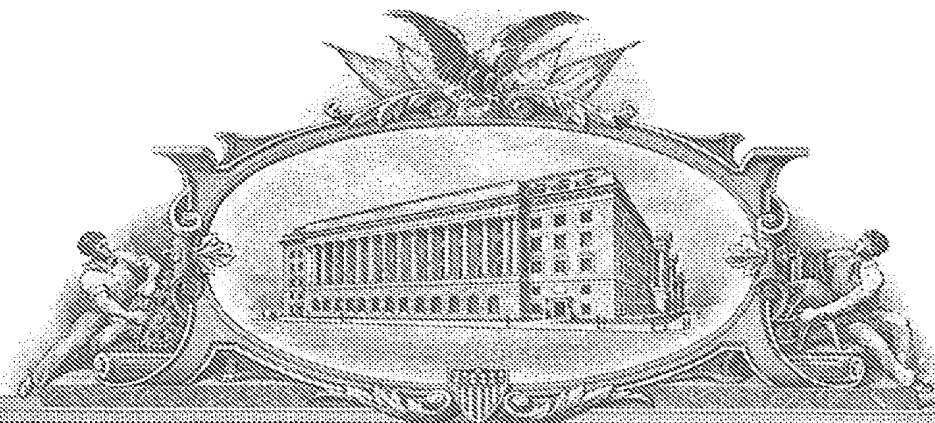

# THE UNITED STATES OF AMERICA

TO ALL TO WHOM THESE PRESENTS SHALL COME:

UNITED STATES DEPARTMENT OF COMMERCE

United States Patent and Trademark Office

*May 02, 2008*

**THIS IS TO CERTIFY THAT ANNEXED HERETO IS A TRUE COPY FROM THE RECORDS OF THE UNITED STATES PATENT AND TRADEMARK OFFICE OF THOSE PAPERS OF THE BELOW IDENTIFIED PATENT APPLICATION THAT MET THE REQUIREMENTS TO BE GRANTED A FILING DATE.**

**APPLICATION NUMBER: 60/504,774**

**FILING DATE: *September 22, 2003***

**RELATED PCT APPLICATION NUMBER: *PCT/US04/30977***

**THE COUNTRY CODE AND NUMBER OF YOUR PRIORITY APPLICATION, TO BE USED FOR FILING ABROAD UNDER THE PARIS CONVENTION, IS *US60/504,774***



Certified by

Under Secretary of Commerce  
for Intellectual Property  
and Director of the United States  
Patent and Trademark Office

02527 U.S. PTO  
09/22/03

4-74

PTO/SB/16 (08-03)  
Approved for use through 07/31/2006. OMB 0651-0032  
U.S. Patent and Trademark Office; U.S. DEPARTMENT OF COMMERCE  
Under the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it displays a valid OMB control number.

**PROVISIONAL APPLICATION FOR PATENT COVER SHEET**

This is a request for filing

37 CFR 1.53(c)

Express Mail Label No.

EV270764558US

22151 U.S. PTO  
60/504774

09/22/03

INVENTOR(S)		
Given Name (first and middle (if any))	Family Name or Surname	Residence (City and either State or Foreign Country)
James E. Stephen W.	Polli Hoag	Ellicott City, MD Baltimore, MD
Additional inventors are being named on the _____ separately numbered sheets attached hereto		
TITLE OF THE INVENTION (500 characters max)		
Method of Assuring Drug and Drug product Identity		
Direct all correspondence to: CORRESPONDENCE ADDRESS		
<input checked="" type="checkbox"/> Customer Number:	33758	
OR		
<input type="checkbox"/> Firm or Individual Name		
Address		
Address		
City	State	Zip
Country	Telephone	Fax
ENCLOSED APPLICATION PARTS (check all that apply)		
<input checked="" type="checkbox"/> Specification Number of Pages	18	<input type="checkbox"/> CD(s), Number
<input type="checkbox"/> Drawing(s) Number of Sheets		<input checked="" type="checkbox"/> Other (specify) Return Postcard
<input type="checkbox"/> Application Date Sheet. See 37 CFR 1.76		
METHOD OF PAYMENT OF FILING FEES FOR THIS PROVISIONAL APPLICATION FOR PATENT		
<input checked="" type="checkbox"/> Applicant claims small entity status. See 37 CFR 1.27.	FILING FEE Amount (\$)	
<input type="checkbox"/> A check or money order is enclosed to cover the filing fees.	80.00	
<input checked="" type="checkbox"/> The Director is hereby authorized to charge filing fees or credit any overpayment to Deposit Account Number: 21-0684		
<input type="checkbox"/> Payment by credit card. Form PTO-2038 is attached.		
The invention was made by an agency of the United States Government or under a contract with an agency of the United States Government.		
<input checked="" type="checkbox"/> No.		
<input type="checkbox"/> Yes, the name of the U.S. Government agency and the Government contract number are:		

[Page 1 of 2]

Respectfully submitted,

SIGNATURE

TYPED or PRINTED NAME AARON D. ADAMS

TELEPHONE 410-706-1933

Date 9/18/03

REGISTRATION NO. 50,278

(if appropriate)

Docket Number: JP-2004-005

**USE ONLY FOR FILING A PROVISIONAL APPLICATION FOR PATENT**

This collection of information is required by 37 CFR 1.51. The information is required to obtain or retain a benefit by the public which is to file (and by the USPTO to process) an application. Confidentiality is governed by 35 U.S.C. 122 and 37 CFR 1.14. This collection is estimated to take 8 hours to complete, including gathering, preparing, and submitting the completed application form to the USPTO. Time will vary depending upon the individual case. Any comments on the amount of time you require to complete this form and/or suggestions for reducing this burden, should be sent to the Chief Information Officer, U.S. Patent and Trademark Office, U.S. Department of Commerce, P.O. Box 1450, Alexandria, VA 22313-1450. DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADDRESS. SEND TO: Mail Stop Provisional Application, Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450.

If you need assistance in completing the form, call 1-800-PTO-9199 and select option 2.

This PTO did not receive the following listed items (s) Small Entity Status

## Method of Assuring Drug and Drug Product Identity

### Field of the Invention

This invention relates generally to a method for assuring drug and drug product identity as a product makes its way from the manufacturer to the retailer.

Within the distribution system of drugs and drug products, there is a need to assure drug and drug product identity. Example scenarios include the need to identify and differentiate authentic and counterfeit drug products, and a need to assure the dispensing of the correct drug product to patients from pharmacies. This invention entails methods to monitor drug product identity of drug products within the drug product distribution system, and achieves this goal through one or more technologies. These technologies and approaches make for a rapid and accurate approach to assure drug and drug product identity, as well as a method to evade drug product from being counterfeited.

Our approach is to apply NIR spectroscopy with a manufacturing method that provides for a dynamic tag system. Future tags are not anticipatable, and perhaps not practically detectable since formulation component(s) themselves provide the tag. The manufacturing method makes use of regulatory mechanisms that were implemented in 1995 (immediate release products) and 1997 (modified release products), in order for the tag system to be dynamic, yet viable from a regulatory point of view.

NIR spectrometry is advantageous in terms of time and disposables. Analysis times are very short (e.g. 1 sec). It is sensitive to multi-component variables, as planned in the described approach. There is essentially no sample preparation. NIR is noninvasive and nondestructive. No reagents are required. Detection limits can be very low.

### Description of the Invention

Counterfeit drug products are an increasing problem, particularly since counterfeiting technology has improved. There remains a need to detect counterfeit drugs, including counterfeit drugs that are similar to the counterfeited products. The invention will include but not limited to the use of near-infrared (NIR), Raman and LIF (laser induced fluorescence) spectroscopy (all possible methods will hence be referred to as spectral methods or NIR methods) to the final dosage form and/or packaging (1-36). Advantages of NIR spectroscopy include its non-invasiveness, potential for low detection limits, and rapidity of analysis (approximately seconds), including minimal or no sample preparation. The vast majority of components in pharmaceutical products exhibit a NIR spectrum.

Either or both the pharmaceutical packaging or dosage form can serve as taggants. For example, embossing, imprinting, printing, coating, dosage form size, and other identification methods can be applied to change the physical appearance of the dosage form (e.g. subtle changes in logo, use of an ultraviolet-dependent dye). Such approach can be used alone or in combination with other modifications of packaging and/or changes in the dosage form (37, 38).

In particular, the components in the dosage form can be varied, approximately batch-to-batch or with some other frequency, to provide a spectral signature for the product. Components of an oral solid dosage forms (e.g. tablets or capsules) include the drug, drug impurities, drug degradants, fillers, disintegrants, binders, lubricants, glidants, colorants, flavoring agents, and coating materials. Some or all of the component levels can be modified, approximately batch-to-batch or with some other frequency, to yield a NIR spectra for the batch or set of batches. The NIR spectra would not be identical for all product batches. A certain batch of product, or certain set of product batches, will have a unique composition, and hence a unique NIR spectra. Only the manufacturer (or agents of the manufacturer) will know the composition and NIR spectra of products from a particular batch. A suspect product can be subjected to NIR analysis and cross-referenced against the authentic NIR spectra. The association between authentic product's batch number and its NIR spectra, along with the ease of measuring NIR spectra, provides a basis to combat counterfeit drugs. Of added benefit, the association between NIR spectra and batch number(s) need not be provided to regulatory agencies or enforcement officials or health care providers, who could still perform field sampling and relay NIR spectra to the manufacturer. Furthermore, past compositions (i.e. past NIR spectra) of previous batches would not be indicative of future compositions (i.e. future NIR spectra). Hence, a counterfeit effort would have no target product to counterfeit, without detection.

This approach employs the NIR region of the electromagnetic spectrum. The NIR region typically includes wavelengths between 700 nm (near the red in the visible spectrum) and 3000 nm (near the infrared stretches of organic compounds). NIR absorbance peaks originate from overtones and combinations of the fundamental (mid-IR) bands and from electronic transitions in the atoms. C-H, N-H, and O-H bonds are responsible for most of the major absorbances. NIR spectrometry is used chiefly to identify or quantify molecules, including unique hydrogen atoms. NIR spectrometry is used to analyze for water, alcohols, amines, and any compounds containing C-H, N-H, and/or O-H groups. Many other bond combinations also provide NIR absorbance peaks

While a pharmaceutical product typically has only one formula, and while manufacturers typically avoid manufacturing changes, this new approach to combat counterfeiting relies on the availability of several (i.e. more than one) formulas for the marketed product. Fortunately, the Food and Drug Administration (FDA) allows for a range of component and composition changes in the manufacturing of products, without onerous regulatory requirements. In the case of immediate release and modified release oral solid dosage forms, changes are denoted Level 1, Level 2, and Level 3 type changes (1,2). Tables 1-9 describe these manufacturing changes. Level 1 changes are those that are unlikely to have any detectable impact on formulation quality and performance; regulatory filing documentation of a Level 1 change is limited to the Annual Report. Level 2 changes are those that could have a significant impact of formulation quality and performance. Level 3 changes are those that are likely to have a significant impact of formulation quality and performance. Tests and filing documentation for a Level 2 change and a Level 3 change each vary depending upon three factors: drug therapeutic range, drug solubility, and drug permeability. Level 1 changes require less burdensome regulatory filing documentation, and represent an advantageous approach than Levels 2 or 3 to vary formulation. Hence, Level 1 changes are a preferable approach to tag authentic product, in order to evade

counterfeiting and facilitate the detection of counterfeiting through NIR spectroscopy. This approach avoids the use of a taggant that is fixed, or is one which is included in the formulation for the sole purpose as a taggant. Our approach to use the formulation's components themselves facilitates the tagging effort, and does so in a more subtle fashion, such that this tagging effort is less detectable and hence counterfeitable.

For example, for each immediate and modified release oral solid dosage forms, a filler can be modified by as much as 5% to provide a NIR spectra that tags the authentic product, and still qualifies as a Level 1 change. Given the ability of NIR spectroscopy to resolve 1% and smaller differences in formulation, a unique tag can be fabricated by varying the filler's level. Moreover, the number of unique NIR signatures can be generated in a multiplicative fashion by modulating two or more components (e.g. vary filler, disintegrant, and binder). Varying filler over 11 levels, disintegrant over 7 levels, and binder over 3 level can results in 231 unique NIR spectra, or more.

While NIR spectroscopy has some previous limited application in pharmaceutical analysis, it does not appear to be applied to combat counterfeit drugs. Our approach employs the formulation itself to provide a dynamic tag system and NIR spectroscopy. Our approach does not employ a fixed tag or a tag whose sole function is to serve as a tag, and thus subject to counterfeiting.

Our approach to apply NIR spectroscopy with a manufacturing method that provides for a dynamic tag system is novel. Future tags are not anticipatable, and perhaps not practically detectable since formulation component(s) themselves provide the tag. The manufacturing method makes use of regulatory mechanisms that were implement in 1995 (immediate release products) and 1997 (modified release products), in order for the tag system to be dynamic, yet viable from a regulatory point of view.

NR spectrometry is advantageous in terms of time and disposables. Analysis times are very short (e.g. 1 sec). It is sensitive to multi-component variables, as planned in the described approach. There is essentially no sample preparation. NIR is noninvasive and nondestructive. No reagents are required. Detection limits can be very low.

The following example is presented in order to more fully illustrate the preferred embodiments of the invention. It should in no way be construed, however, as limiting the broad scope of the invention.

#### Example 1

**Materials.** The following drug substances and excipients were used as received: aspirin (Spectrum, Gardena; CA), prednisone (Sigma; St Louis, MO), indomethacin (Spectrum; Gardina, CA), acyclovir (Spectrum; Gardena, CA), microcrystalline cellulose (Emocel 90M, Mendell; Patterson, NY), magnesium stearate (Spectrum, Gardina, CA), croscarmellose sodium (FMC Biopolymer; Princeton, NJ), starch (Lycatab C, Roquette; Lestrem, France), and lactose monohydrate (Super-tab, The Lactose Company; Hawera, New Zealand),

**Formulation Methods.** Three tablet formulations were designed and evaluated for each of four drugs, such that 12 formulations were made. The four drugs were aspirin, prednisolone, indomethacin, and acyclovir, and are denoted as drug A, B, C, and D, respectively. The drugs differ in their therapeutic uses, physicochemical properties, spectral properties, and dose ranges. For each drug, three tablet formulations were fabricated. Tables 1-4 describe the composition of the 12 formulations and refer to formulations A1, A2, A3, B1, etc. In each table, the first formulation is denoted the reference formulation (i.e. A1, B1, C1, and D1 are reference formulations). For each drug, the formulations were varied within the SUPAC IR level 1 tolerances by varying one or more excipients, relative to the reference formulation, resulting in the second and third formulations (i.e. formulations A2 and A3 were variants for formulation A1; formulations B2 and B3 were variants for formulation B1).

Variant formulations were attained through the following changes, relative to the reference. For aspirin, microcrystalline cellulose was increased and decreased. For prednisone, magnesium stearate was increased and decreased. For indomethacin, microcrystalline cellulose and croscarmellose sodium were simultaneously varied. For acyclovir, microcrystalline cellulose and lactose monohydrate were simultaneously varied. In some cases the tablet weight changed.

**Near-IR Methods.** The formulations were scanned and analyzed by Foss NIRSystems Rapid Content Analyzer™. The following test conditions were used. Samples were placed into sealed glass scintillation vials and scanned in reflectance mode; each sample was scanned 62 times and averaged into one spectrum; the wavelength range was 400 nm to 2500 nm with samples collected every 2 nm. The raw spectral data were converted into absorbance and 2<sup>nd</sup> derivative values using Foss's Vision software package.

Also, because this device can detect the drug type and its dose this could also be used to greatly minimize the possibility of the pharmacist dispensing the wrong drug and the wrong dose of a drug

A unique NIR signature could be engineered into the packaging, much like a certificate of authenticity on a CD or commercial software package. A company could have several hundred types and assign each one with a lot number.

It should also be noted that varying the formulation to evade and detect counterfeiting is one approach in the application of spectral methods. Another approach does not make use of the intentional periodic variation of the formulation as described above.

Various applications for this technology include, but should not be limited to

- 1) methods to evade and detect counterfeit drug products (and counterfeit drug substances and counterfeit excipients)
- 2) methods to assure drug product distribution integrity (and drug substance integrity and excipient integrity)
  - a. levels for monitoring drug product distribution:
    - a1. pharmaceutical manufacturers
    - a2. pharmaceutical wholesalers
    - a3. pharmaceutical distributors

a4. pharmacies (NOTE: While others may be primarily interested in detecting counterfeit drug products, pharmacies would also be interesting in assuring the correct product is being dispensed.)

a5. pharmaceutical repackagers

a6. FDA field monitoring, as well as regulators in other countries

Typical example applications include but are not limited to

a. Evade and detect counterfeit drug products by FDA field inspectors and/or health care worker (e.g. pharmacist, nurse) working with manufacturer of the authentic product; FDA field inspectors and/or health care workers obtain NIR spectrum of suspect products and relay the spectrum data to the manufacturer of the authentic product. Agents of the manufacturer of the authentic product may also inspect samples in the field by obtaining NIR spectrum of suspect products.

b. A pharmacy employs NIR to assure that the dispensing robot dispenses the correct product.

c. A pharmacy employs NIR to assure the dispensed product is the correct product (e.g. NIR built into tablet/capsule counter, a semi-automated dispensing device).  
(e.g. to avoid accidental overdose or misadventure in pharmacy dispensing)

#### Brief Description of the Figures

Table 1. Level 1 Component and Composition Changes for Immediate Release Oral Solid Dosage Forms

Table 2. Level 2 Component and Composition Changes for Immediate Release Oral Solid Dosage Forms

Table 3. Level 3 Component and Composition Changes for Immediate Release Oral Solid Dosage Forms

Table 4. Level 1 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (nonrelease controlling excipient)

Table 5. Level 2 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (nonrelease controlling excipient)

Table 6. Level 3 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (nonrelease controlling excipient)

Table 7. Level 1 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (release controlling excipient)

Table 8. Level 2 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (release controlling excipient)



Table 9. Level 3 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (release controlling excipient)

Table 10. Schematic of possible areas of use for this invention within the commercial pipeline

Table 11. Composition of Aspirin Formulations

Table 12. Composition of Prednisone Formulations

Table 13. Composition of Indomethacin Formulations

Table 14. Compositions of Acyclovir Formulations

Table 15. 2<sup>nd</sup> Derivative of Absorbance vs. Wavelength: Aspirin Formulations

Table 16. 2<sup>nd</sup> Derivative of Absorbance vs. Wavelength: Prednisone Formulations

Table 17. 2<sup>nd</sup> Derivative of Absorbance vs. Wavelength: Indomethacin Formulations

Table 18. 2<sup>nd</sup> Derivative of Absorbance vs. Wavelength: Acyclovir Formulations

## References

The following publications, which may be cited above, are hereby incorporated by reference in their entireties for all purposes:

1. Rubinovitz, R. 2003. Powder blend uniformity by NIR. *Manufacturing Chemist* 74: 37-38.
2. Beyer, J.; Steffens, K. J. 2003. Calibration models for determination of water content in pharmaceutical excipients using near-infrared spectroscopy (NIRS). *Pharmazeutische Industrie* 65: 186-192.
3. Kuny, T.; Schatz, C.; Ulmschneider, M.; Marrer, S.; Leuenberger, H. 2003. Non-destructive dissolution testing correlation. *Dissolution Technology* 10: 22-24,26,28.
4. Gustafsson, C.; Nystroem, C.; Lennholm, H.; Bonferoni, M. C.; Caramella, C. M. 2003. Characteristics of hydroxypropyl methylcellulose influencing compactibility and prediction of particle and tablet properties by infrared spectroscopy. *Journal of Pharmaceutical Sciences* 92: 494-504.
5. Anonymous 2002. NIR enters the pharma arena. *Manufacturing Chemist* 73: 59-60.
6. Ochs, R. 2002. Sounding the alarm on counterfeit drugs. *Pharmacy Review* 26: 34-35.
7. Yang, H.; Irudayaraj, J. 2002. Rapid determination of vitamin c by NIR, MIR and FT-Raman techniques. *Journal of Pharmacy & Pharmacology* 54: 1247-1255.
8. Jorgensen, A., et al. 2002. Hydrate formation during wet granulation studied by spectroscopic methods and multivariate analysis. *Pharmaceutical Research* 19: 1285-1291.

9. Reich, G. 2002. Potential of attenuated total reflection infrared and near-infrared spectroscopic imaging for quality assurance/quality control of solid pharmaceutical dosage forms. *Pharmazeutische Industrie* 64: 870-874.
10. Scott, P. 2002. Process analytical technology: Applications to the pharmaceutical industry. *Dissolution Technology* 9: 6,8.
11. Rager, I.; Roos, G.; Schmidt, P. C.; A., K. K. 2002. Rapid quantification of constituents in St. John's wort extracts by NIR spectroscopy. *Journal of Pharmaceutical & Biomedical Analysis* 28: 439-446.
12. Cl, J. 2002. Experts give clues on how to stay clear of counterfeit drugs. *Drug Topics* 146: 66.
13. Ritchie, G. E., et al. 2002. Validation of a near-infrared transmission spectroscopic procedure Part B: Application to alternate content uniformity and release assay methods for pharmaceutical solid dosage forms. *Journal of Pharmaceutical & Biomedical Analysis* 29: 159-171.
14. El-Hagrasy, A. S.; Morris, H. R. D.; Amico, F.; Lodder, R. A.; Drennen, J. K. 2001. Near-infrared spectroscopy and imaging for the monitoring of powder blend homogeneity. *Journal of Pharmaceutical Sciences* 90: 1298-1307.
15. Harris, S. C.; Walker, D. S. 2000. Quantitative real time monitoring of dryer effluent using fiber optic near infrared spectroscopy. *Journal of Pharmaceutical Sciences* 89: 1180-1186.
16. Anonymous 2001. NIR spectroscopy in a manufacturing plant. *Manufacturing Chemist* 72: 23-24.
17. Rein, H. 2000. NIR-Vis spectroscopy: modern analysis process. *Deutsche Apotheker Zeitung* 140: 45-50, 57-58.
18. Rantanen, J., et al. 2000. In line moisture measurement during granulation with a four wavelength near infrared sensor: evaluation of particle size and binder effects. *European Journal of Pharmaceutics & Biopharmaceutics* 50: 271-276.
19. Brashear, R. L.; Flanagan, D. R.; Luner, P. E.; Seyer, J. J.; Kemper, M. S. 1999. Diffuse reflectance near infrared spectroscopy as a nondestructive analytical technique for polymer implants. *Journal of Pharmaceutical Sciences* 88: 1348-1353.
20. Candolfi, A.; De Maesschalck, R.; Massart, D. L.; Hailey, P. A.; Harrington, A. C. 1999. Identification of pharmaceutical excipients using NIR spectroscopy and SIMCA. *Journal of Pharmaceutical & Biomedical Analysis* 19: 923-935.
21. Glover, M.; Lewis, T. 1999. Simplifying tablet validation analysis. *Manufacturing Chemist* 70: 70-72.
22. Buckton, G.; Yonemochi, E.; Hammond, J.; Moffat, A. 1998. Use of infrared spectroscopy to detect changes in the form of amorphous and crystalline lactose. *International Journal of Pharmaceutics* 168: 231-241.
23. Gold, T. B.; Buice, R. G.; Lodder, R. A.; Digenis, G. A. 1998. Detection of formaldehyde-induced crosslinking in soft elastic gelatin capsules using near-infrared spectrophotometry. *Pharmaceutical Development & Technology* 3: 209-214.
24. Candolfi, A.; Wu, W.; Massart, D. L.; Heuerding, S. 1998. Comparison of classification approaches applied to NIR-spectra of clinical study lots. *Journal of Pharmaceutical & Biomedical Analysis* 16: 1329-1347.
25. Ren, Y. L.; Li, W. 1998. Analysis of sulfamethoxazole (SZ) by PLS-NIR diffuse reflectance spectrophotometry. *Chinese Journal of Pharmaceutical Analysis* 18: 30-33.

26. Adams, D.; Brown, G. P.; Fritz, C.; Todd, T. R. 1998. Calibration of a near-infrared (NIR) H<sub>2</sub>O<sub>2</sub> vapor monitor. *Pharmaceutical Engineering* 18: 66-68, 70-76, 78-82.
27. Gold, T. B.; Buice, R. G.; Lodder, R. A.; Digenis, G. A. 1997. Determination of extent of formaldehyde induced crosslinking in hard gelatin capsules by near infrared spectrophotometry. *Pharmaceutical Research* 14: 1046-1050.
28. Higgins, M. 1997. Pinpointing production problems with NIR analysis. *Manufacturing Chemist* 68: 38-39.
29. Buchanan, B. R.; Baxter, M. A.; Chen, T. S.; Qin, X. Z.; Robinson, P. A. 1996. Use of near-infrared spectroscopy to evaluate an active in a film coated tablet. *Pharmaceutical Research* 13: 616-621.
30. Buice, R. G.; Gold, T. B.; Lodder, R. A.; Digenis, G. A. 1995. Determination of moisture in intact gelatin capsules by near-infrared spectrophotometry. *Pharmaceutical Research* 12: 161-163.
31. Van Der Vlies, C.; Kaffka, K. J.; Plugge, W. 1995. Qualifying pharmaceutical substances by fingerprinting with NIR spectroscopy and PQS. *Pharmaceutical Technology International* 7: 43-44, 46, 48, 56.
32. Galante, L. J.; Brinkley, M. A.; Lodder, R. A. 1992. Bacterial monitoring in vials using a spectrophotometric assimilation method. *Pharmaceutical Research* 9: 357-364.
33. Zannikos, P. N.; Li, W. I.; Drennen, J. K.; Lodder, R. A. 1991. Spectrophotometric prediction of the dissolution rate of carbamazepine tablets. *Pharmaceutical Research* 8: 974-978.
34. Drennen, J. K.; Lodder, R. A. 1990. Nondestructive near-infrared analysis of intact tablets for determination of degradation products. *Journal of Pharmaceutical Sciences* 79: 622-627.
35. Kamat, M. S.; Lodder, R. A.; DeLuca, P. P. 1989. Near-infrared spectroscopic determination of residual moisture in lyophilized sucrose through intact glass vials. *Pharmaceutical Research* 6: 961-965.
36. Loder, R. E. 1977. Use of hyperbaric oxygen in paralytic ileus. *British Medical Journal* 1: 1448-1449.
37. Anonymous 2003. Faking it. *Manufacturing Chemist* 74: 23-24.
38. Karsten, R. 2003. Tags tag counterfeits. *Manufacturing Chemist* 74: 26-27.

Table 1. Level 1 Component and Composition Changes for Immediate Release Oral Solid Dosage Forms

Excipient	Percent Excipient (w/w) Out of Total Target Dosage Form Weight
Filler	+/- 5%
Disintegrant	+/- 3%
starch	+/- 1%
other	+/- 0.5%
Binder	+/- 0.25%
Lubricant	+/- 1%
calcium or magnesium stearate	
other	
Glidant	+/- 1%
talc	+/- 0.1%
other	+/- 1%
Film coat	

Table 2. Level 2 Component and Composition Changes for Immediate Release Oral Solid Dosage Forms

Excipient	Percent Excipient (w/w) Out of Total Target Dosage Form Weight
Filler	+/- 10%
Disintegrant	+/- 6%
starch	+/- 2%
other	+/- 1%
Binder	+/- 0.5%
Lubricant	+/- 2%
calcium or magnesium stearate	
other	
Glidant	+/- 2%
talc	+/- 0.2%
other	+/- 2%
Film coat	

Table 3. Level 3 Component and Composition Changes for Immediate Release Oral Solid Dosage Forms

Excipient	Percent Excipient (w/w) Out of Total Target Dosage Form Weight
Filler	Greater than +/- 10%
Disintegrant	
starch	Greater than +/- 6%
other	Greater than +/- 2%
Binder	Greater than +/- 1%
Lubricant	
calcium or magnesium stearate	Greater than +/- 0.5%
other	Greater than +/- 2%
Glidant	
talc	Greater than +/- 2%
other	Greater than +/- 0.2%
Film coat	Greater than +/- 2%

Table 4. Level 1 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (nonrelease controlling excipient)

Excipient	Percent Excipient (w/w) Out of Total Target Dosage Form Weight
Filler	+/- 5%
Disintegrant	
starch	+/- 3%
other	+/- 1%
Binder	+/- 0.5%
Lubricant	
calcium or magnesium stearate	+/- 0.25%
other	+/- 1%
Glidant	
talc	+/- 1%
other	+/- 0.1%
Film coat	+/- 1%

Table 5. Level 2 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (nonrelease controlling excipient)

Excipient	Percent Excipient (w/w) Out of Total Target Dosage Form Weight
Filler	+/- 10%
Disintegrant	+/- 6%
starch	+/- 2%
other	+/- 1%
Binder	+/- 1%
Lubricant	+/- 0.5%
calcium or magnesium stearate	+/- 2%
other	+/- 2%
Glidant	+/- 2%
talc	+/- 0.2%
other	+/- 2%
Film coat	+/- 2%

Table 6. Level 3 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (nonrelease controlling excipient)

Excipient	Percent Excipient (w/w) Out of Total Target Dosage Form Weight
Filler	Greater than +/- 10%
Disintegrant	Greater than +/- 6%
starch	Greater than +/- 2%
other	Greater than +/- 1%
Binder	Greater than +/- 1%
Lubricant	Greater than +/- 0.5%
calcium or magnesium stearate	Greater than +/- 2%
other	Greater than +/- 2%
Glidant	Greater than +/- 2%
talc	Greater than +/- 0.2%
other	Greater than +/- 2%
Film coat	Greater than +/- 2%

Table 7. Level 1 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (release controlling excipient)

Excipient	Percent Excipient (w/w) Out of Total Release Controlling Excipient Content in the Modified Release Solid Oral Dosage Form
Any release controlling excipient(s)	+/- 5%

Table 8. Level 2 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (release controlling excipient)

Excipient	Percent Excipient (w/w) Out of Total Release Controlling Excipient Content in the Modified Release Solid Oral Dosage Form
Any release controlling excipient(s)	+/- 10%

Table 9. Level 3 Component and Composition Changes for Modified Release Oral Solid Dosage Forms (release controlling excipient)

Excipient	Percent Excipient (w/w) Out of Total Release Controlling Excipient Content in the Modified Release Solid Oral Dosage Form
Any release controlling excipient(s)	Greater than +/- 10%

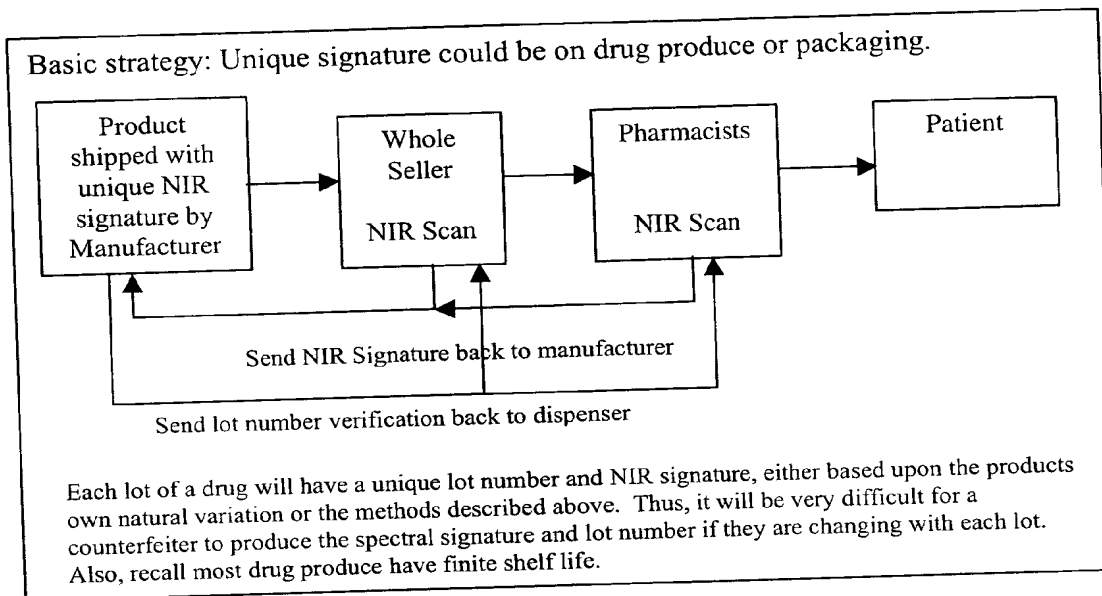




Table 1. Composition of Aspirin Formulations

Component	Formulation A1 (mg/tab)	Formulation A2 (mg/tab)	Formulation A3 (mg/tab)
Aspirin	325	325	325
Microcrystalline cellulose	73	83	63
Magnesium stearate	2	2	2
TOTAL WEIGHT	400	410	390

Table 2. Composition of Prednisone Formulations

Component	Formulation B1 (mg/tab)	Formulation B2 (mg/tab)	Formulation B3 (mg/tab)
Prednisone	5	5	5
Microcrystalline cellulose	94.5	94.5	94.5
Magnesium stearate	0.5	0.75	0.25
TOTAL WEIGHT	100	100.25	99.75


Table 3. Composition of Indomethacin Formulations

Component	Formulation C1 (mg/tab)	Formulation C2 (mg/tab)	Formulation C3 (mg/tab)
Indomethacin	25	25	25
Microcrystalline cellulose	71.5	74	69
Croscarmellose sodium	3	2	4
Magnesium stearate	0.5	0.5	0.5
TOTAL WEIGHT	100	101.5	98.5

Table 4. Composition of Acyclovir Formulations

Component	Formulation D1 (mg/tab)	Formulation D2 (mg/tab)	Formulation D3 (mg/tab)
Acyclovir	200	200	200
Microcrystalline cellulose	113.26	120.26	106.26
Starch	35	27.99	41.99
Magnesium stearate	1.75	1.75	1.75
TOTAL WEIGHT	350	350	350


## 2<sup>nd</sup> Derivative of Absorbance Versus Wavelength: Aspirin Formulations

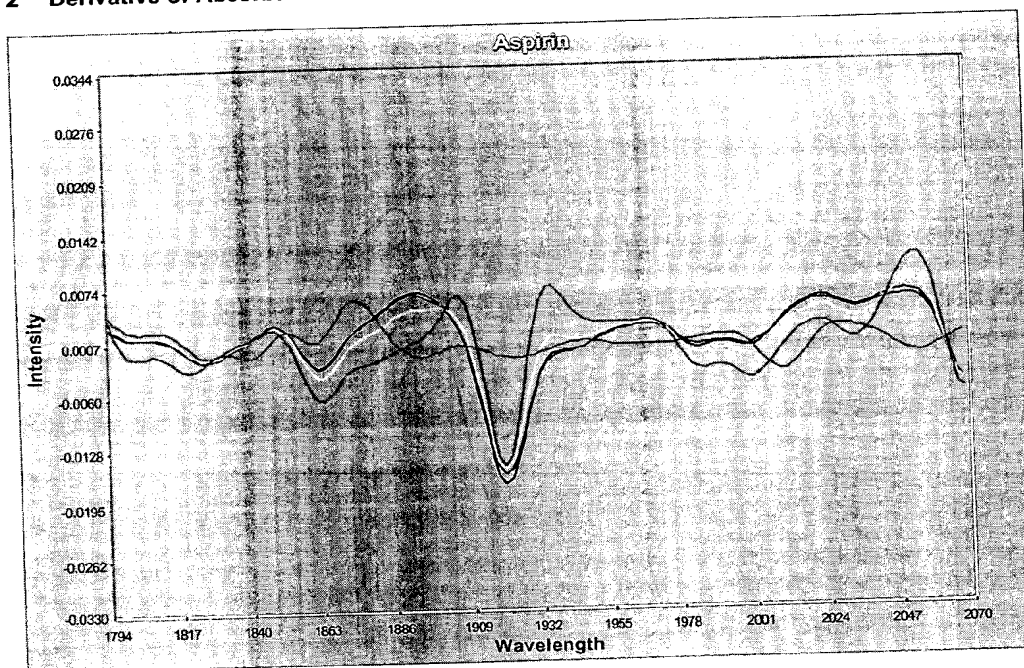

2<sup>nd</sup> Derivative of Absorbance Versus Wavelength: Aspirin Formulations. Formulations A3 (Yellow), A1 (Blue), and A2 (Red) contained increasing amounts of microcrystalline cellulose. The intensities around 1995 nm and 2055 nm reflect NIR to differentiate the formulations. The profiles of pure microcrystalline cellulose (Light Blue) and pure aspirin (Green) are also shown.

Table 15


**2<sup>nd</sup> Derivative of Absorbance Versus Wavelength: Prednisone Formulations.**

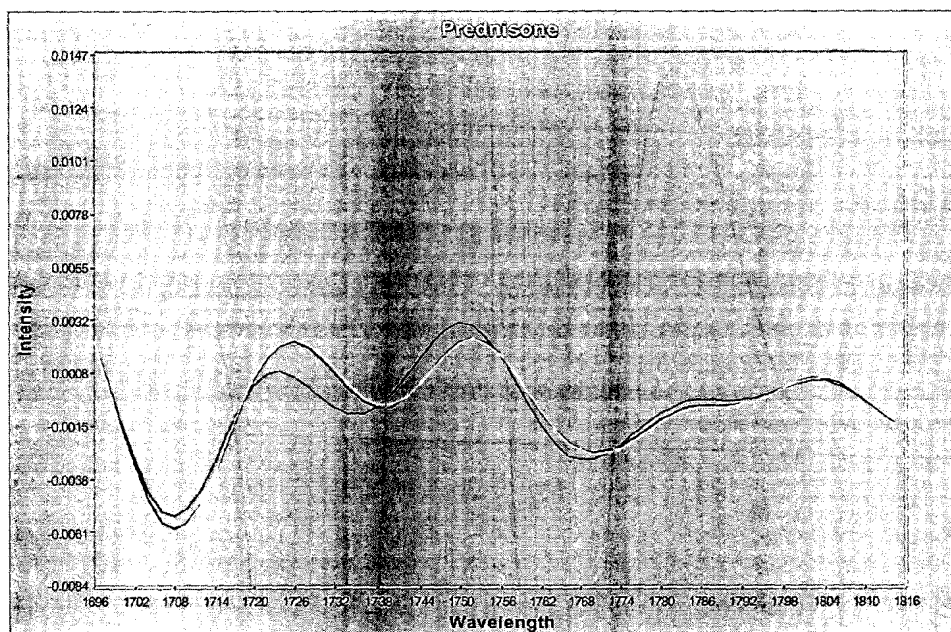

2<sup>nd</sup> Derivative of Absorbance Versus Wavelength: Prednisone Formulations. Formulations B3 (Yellow), B1 (Blue), and B2 (Red) contained increasing amounts of magnesium stearate. The intensities around 1705 nm, as well as the regions between 1725-1735 nm and 1775-1790 nm, reflect NIR to differentiate the formulations. The profile of pure magnesium stearate (Light Blue) is also shown.

Table 16


## 2<sup>nd</sup> Derivative of Absorbance Versus Wavelength: Indomethacin Formulations

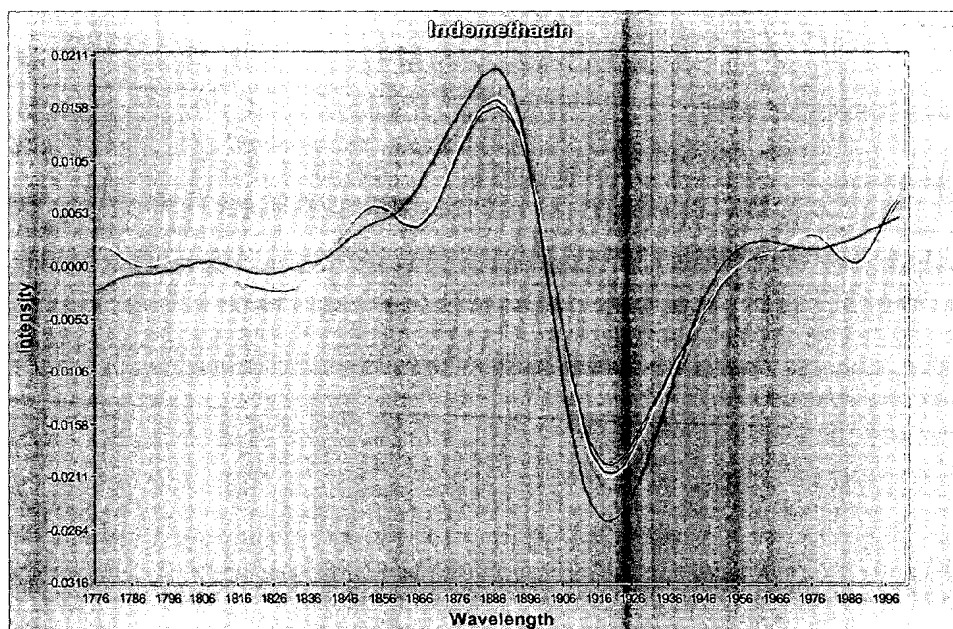

2<sup>nd</sup> Derivative of Absorbance Versus Wavelength: Indomethacin Formulations. Formulations C3 (Yellow), C1 (Blue), and C2 (Red) contained increasing amounts of microcrystalline cellulose, as well as decreasing amounts of croscarmellose sodium. The intensities around 1890 nm and 1920 nm reflect NIR to differentiate the formulations. The profiles of pure microcrystalline cellulose (Light Blue) and pure croscarmellose sodium (Purple) are also shown.

Table 17

**2<sup>nd</sup> Derivative of Absorbance Versus Wavelength: Acyclovir Formulations**

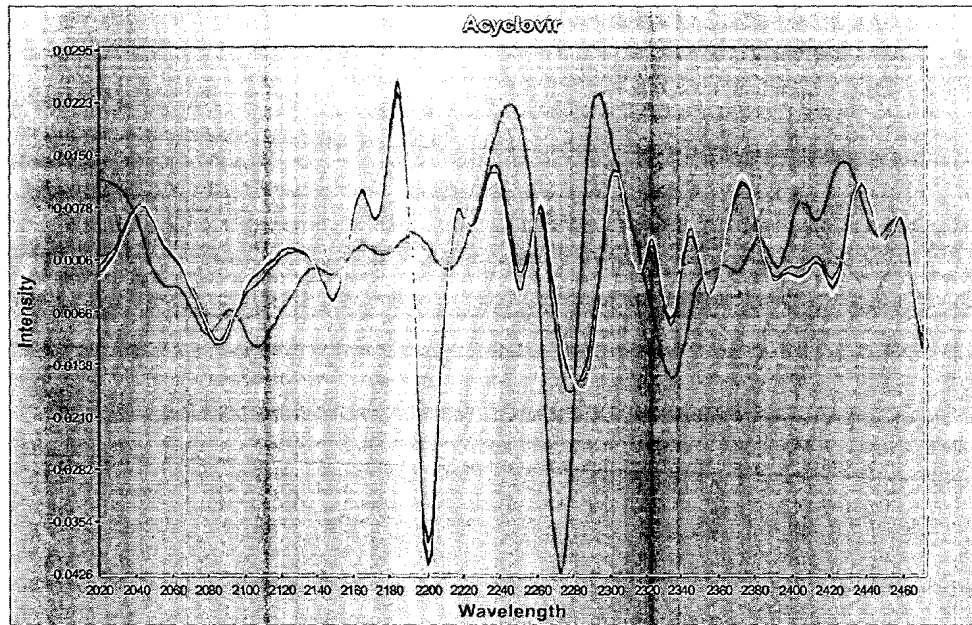

2<sup>nd</sup> Derivative of Absorbance Versus Wavelength: Acyclovir Formulations. Formulations C3 (Yellow), C1 (Blue), and C2 (Red) contained increasing amounts of microcrystalline cellulose, as well as decreasing amounts of starch. The intensities around 2175 nm, 2205 nm, 2225 nm, 2250 nm, 2265 nm, 2320 nm, 2345 nm, 2365 nm, as well as the regions between 2100-2130 nm and 2380-2420 nm, reflect NIR to differentiate the formulations. The profiles of pure microcrystalline cellulose (Gray) and pure starch (Light Blue) are also shown.

Table 18